

IBEX

User Guide

IBEX 2-Step LED Slit Lamp - Compact



Introduction

Congratulations on your new purchase and welcome to IBEX! We appreciate your patronage and the opportunity to exceed your expectations with outstanding products and services. As a valued customer, your experience with IBEX is our top priority, and we strive for a 100% customer satisfaction rate.

Your new IBEX 2-Step LED Slit Lamp - Compact arrives with a comprehensive 3-Year Warranty, should you have any questions or concerns, please do not hesitate to contact us.

Product Registration

If you purchased your product directly through IBEX, **your product is already registered.**

If you purchased your product from one of our many Business Partners, please register your products.

As part of our ongoing efforts to provide you with the best possible customer support and user experience, we urge you to register your new product at <https://www.ibexeye.com/welcome>. *It's fast & easy.*

IBEX LED Slit Lamps are serialized in order to track product versions and revision changes. The serial number may be required to complete any necessary services and or repairs. Registering your product serial number simply allows our IBEX support team to better understand the product revisions in order to expedite your service.



Before you begin using your new **IBEX 2-Step LED Slit Lamp - Compact**, please read, understand and follow all safety instructions. This will insure safe operation, optimum performance and a longer service life for your instrument. Be sure to retain this manual for future reference.

Warning and Information Symbols



General Warning



General Mandatory Action



Hot Surface



Read Operator's Manual



Voltage Hazard



Type B Applied Part



Optical Radiation



Notes



Table of Contents

Introduction	1	Biomicroscope	14
Product Registration	1	Preparing the patient.....	15
Warning and Information Symbols	1	Patient Alignment	15
Table of Contents	2	Basic Instrument Operation.....	16
Product Overview	3	Maintenance & Cleaning	16
Indications for use.....	3	Routine maintenance	16
Major Components	3	Fuse Replacement	17
Illumination Tower and Biomicroscope Swing Arm Assembly	4	General cleaning	17
Slit Lamp Transitional Base Assembly, Table base with headrest	5	Recommended Lens Cleaning Technique and Procedure	18
Accessories	6	Technical Specifications.....	19
Optional Accessories:.....	6	General Classifications.....	19
Installation	7	Basic Troubleshooting	19
Unpacking	7	Block diagram	20
Ophthalmic Stand Installation	8	Storage, Transportation & Environment	20
Leveling the slit lamp	10	Warranty.....	20
Motorized Table Installation	11	Quality Systems	21
Installing Optional Applanation Tonometer	12	Safety Instructions	21
Usage - Getting Started.....	12	Safety Standards	22
Knowing your Slit Lamp.....	13	Standards Compliance	23
Telescoping eye cup adjustment	13	Contact Information	24
Eyepiece setting	13		
PD setting	14		
Slit width control.....	14		
Slit angle control	14		
Setting aperture size	14		
Filter Selection	14		


Product Overview

Indications for use


All **IBEX Slit Lamps** are specialized LED illuminated biomicroscopes that provide a stereoscopic magnified view of the eye structures in crystal clear detail, enabling anatomical diagnoses to be made for a variety of eye conditions. They are precision instruments intended for use, by professionally trained personnel, in the examination of the anterior eye segment from the cornea epithelium to the posterior capsule. They are used to aid in the diagnosis of diseases or trauma in which affects the structural properties of the anterior eye segment.

Major Components

IBEX Slit Lamps consists of 3 major components; **main body**, **headrest assembly** and **table base**.

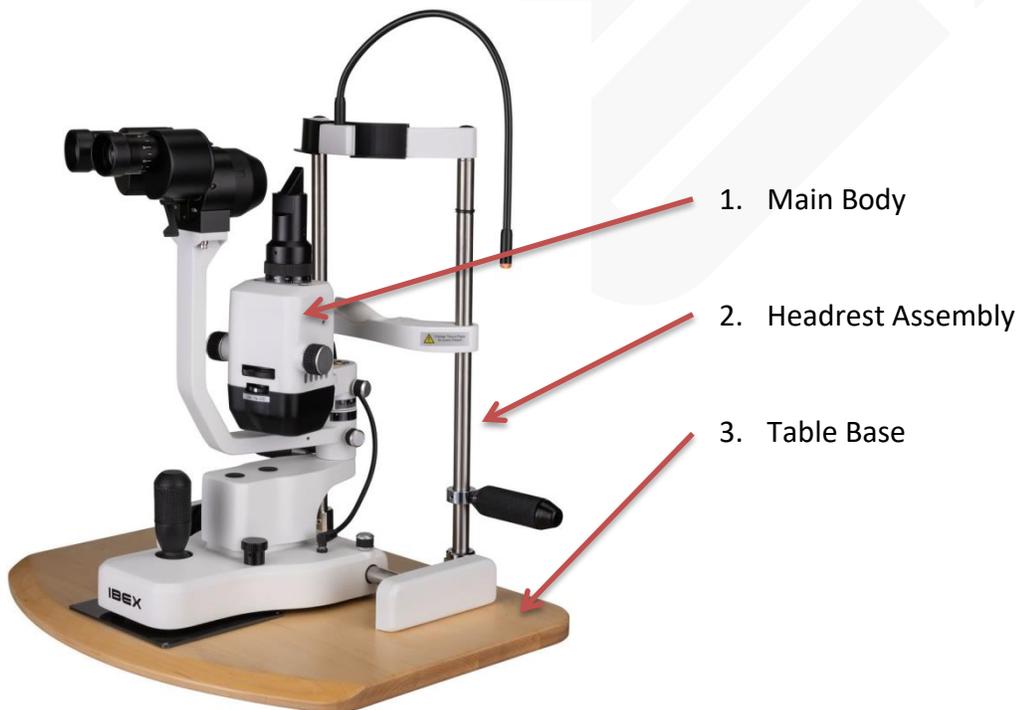

Illumination Tower and Biomicroscope Swing Arm Assembly

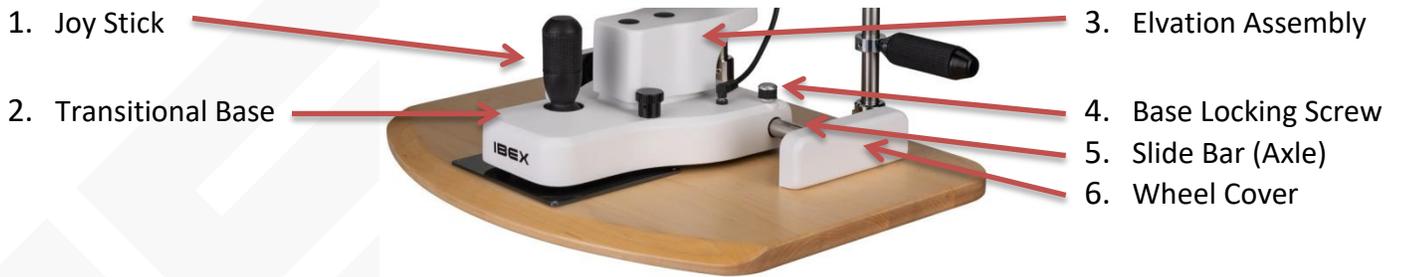

Slit Lamp Transitional Base Assembly



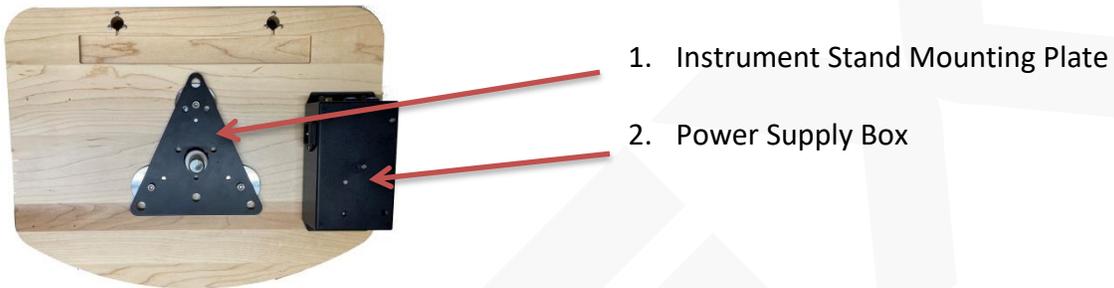


Table base with Headrest

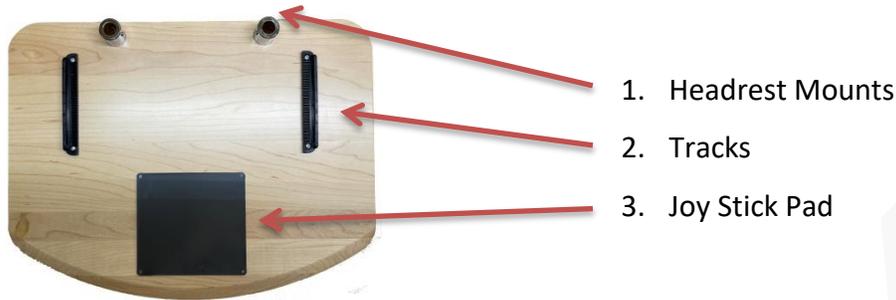

Accessories

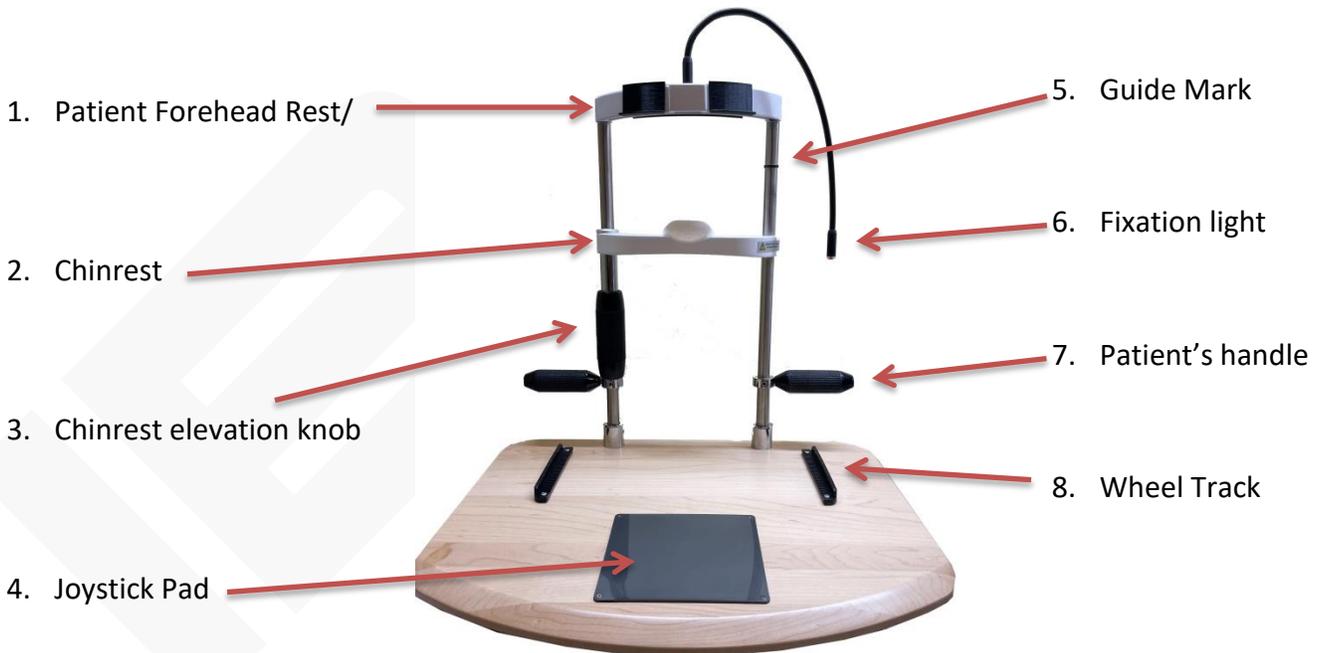

1. Test Rod
2. Chin paper Pins
3. Pivot Cover
4. Track Covers
5. Chin Paper
6. Link Cable

7. Patient's Handles
8. Forehead Rest/Doctor handrest
9. Dust Cover
10. Breath Shield
11. Fuses
12. Allen Wrench Set

Notes: Eyepieces (15X) only come with the 2L model

Optional Accessories: IBEX Model R Tonometer

Installation

Unpacking

Your New **IBEX Slit Lamp** arrived in 1 box, the Slit lamp body, headrest assembly and table top assembly. Slit lamp body is fully assembly for easy installation. They are double boxed and layered in foam for added safety.




Your **IBEX** Slit Lamp can be mounted on a 3rd party instrument stand or a motorized table.



Ensure that your instrument stand or motorized table is on a **level and stable** surface before installation of your new slit lamp. Use the level from the accessory kit if necessary.

Ophthalmic Stand Installation

Step 1: Installation of the table base on your instrument stand.

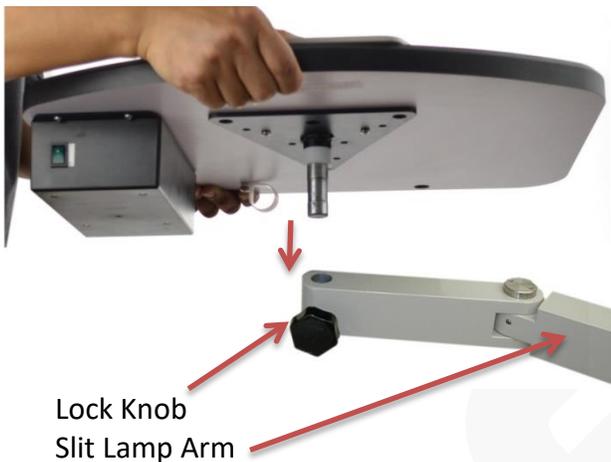

- Remove the table base from shipping carton.
- Insert the mounting post (on the underside of the table base) into the slit lamp arm on the instrument stand.



You may have to loosen the lock knob on the slit lamp arm.

Step 2: Installation of the headrest on the table base.



- Remove headrest from Slit Lamp box.
- Insert the headrest cables through the headrest mounts.
- Install the headrest onto the table base and tighten the four set screws to secure the headrest.


Step 3: Installation of the slit lamp on the table base.



a. Install the slit lamp onto the table base.



b. Insure the slit lamp is aligned evenly on the tracks.

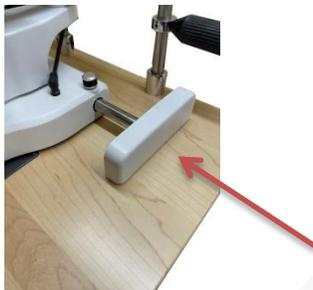

c. Install the track covers by inserting them under the tracks.

Step 4: Connection of cables.



a. Remove link cable from the accessory box and connect to the slit lamp transitional base.

b. Connect the fixation light cable and intensity cable to the back of the power supply.

c. Connect the main power cord into the side of the power supply and to your AC power source.



Leveling the slit lamp

Your slit lamp is come with a smart leveling mechanism. It allows you to level the slit lamp on an instrument stand.



- Step 1: Place the bubble level on the table top in front of the slit lamp gliding base to check the level on x-axis. Then put the bubble level beside the joystick pad to check the level of y-axis.



- Step 2: Loosen the allen screw until the bubble sets in the center of the level. You may need to adjust all 3 allen screws multiple times to achieve centering.



You can make an adjustment of as much as 5° with respect to the table top.

- Step 3: Turn the thumb nut until it against the mounting plate to lock the position.

Motorized Table Installation

Installation of the table base on an **IBEX** motorized table.



Again insure that the surface where you are placing your table is level and stable. Use the level from the accessory kit.

Step 1:

	<ol style="list-style-type: none"> a. Remove table base from the small shipping carton. b. Remove the triangular mounting bracket from the bottom of the table base if applicable.
--	--

Step 2:

	<ol style="list-style-type: none"> c. Place the motorized table mounting plate (if the motorized table came with one) or the motorized table column in the center of the table top. And the controller near the front end of the table top. Mark and drill the mounting holes.
--	---

Step 3:

	<ol style="list-style-type: none"> d. Attach the motorized table column and the up/down controller with the provided screws and washers. Then connect the cable.
--	---

Installing Optional Applanation Tonometer



Step 1: Remove the silver screw from top of the slit lamp biomicroscope head.



Step 2: Attach the Tonometer mounting block onto the biomicroscope head using the countersunk screw that is supplied with the Tonometer.



Step 3: Slide the Tonometer Swing Arm Mount down onto the vertical post of the Mounting Block.



Step 4: Insert the Prism completely into the Prism Holder so that the axis marks on the Prism and Prism Holder meet.

Usage - Getting Started



Before use, the slit lamp should be allowed to adjust to the ambient room temperature for several hours. This is especially important when the unit has been stored or transported in a cold environment.

Before setting up the Biomicroscope eye pieces, set the magnification to 1.6x by flipping the lens magnification changing lever from left to right.

Each eyepiece and the PD must be individually adjusted for each examiner.

Telescoping eye cups at the front of the oculars enable users with or without glasses to see the maximum field of view.

Insert test rod in place of pivot cover.

Turn the test rod so that the flat surface is facing toward the Bio-microscope.

Return both the illumination tower and Bio-microscope swing arms to the central position.



Never attempt to connect or disconnect this instrument with **wet or damp** hands.

Do not disconnect this instrument by pulling on the power cord. Grip and hold the connector firmly when disconnecting from the power source.



LED fixtures can reach high temperatures during use – allow ample time to cool before handling any potentially hot surface.



Switch the slit lamp ON using the rocker switch on the front of the power supply. The rocker switch will light up green when the switch is in the ON position. Slowly turn the intensity control knob on the transitional base until you see the light projected onto the test rod. If you do not see any light on the test rod turn the slit width control knob to full width.

Knowing your Slit Lamp

Telescoping eye cup adjustment

The dual purpose of the eye cups are excluding extraneous sidelight and positioning the eye pupils to the correct distance from the eyepieces so as to see the instrument's full field of view. As general rule of thumb, pull the eye cups out for non-eyeglass wearers and push the eye cups in for eyeglass wearers.

Eyepiece setting

Turn both eyepieces to maximum (+) diopter position. Close your **non-dominant eye** and with your **dominant eye** look through the ocular and focusing on the test rod circle of light. Slowly turn the eyepiece toward the (-) diopter position until the circle of light surface of the test rod is in focus. Repeat the same procedure for the other eyepiece.

PD setting

Begin by holding the binocular head in a normal viewing position. Then, while grasping the ocular barrels firmly with your hands, move them either closer together or spread them farther apart until the images seen by your eyes form a single circular field of view.



Slit width control

Slit width can be changed continuously (from 0 to 14mm depending on the aperture selected) by rotating either of the two knurled slit width control knobs to the desired width.

Slit angle control

By turning the dial clockwise or counterclockwise, the slit image will revolve at any angle in the vertical or horizontal direction. The angle of image can be rotated from 0° and 180° with indicated in 5° increments on the rotation angle scale.

Setting aperture size

Turn the aperture control knob to the desired aperture. The following separate diameters are available: 14, 10, 5, 1, 0.3mm. Width and length of the slit image can be changed continuously from 12mm to 1mm, which is indicated on the aperture dial.

Filter Selection

Filter selection can be changed by turning the filter selection dial clockwise or counterclockwise. The order for the filter in clockwise direction is: Open, Yellow, Cobalt blue, and Green as indicated by the colored dot on the filter dial.

Biomicroscope

The magnification can be changed from 1x to 1.6x, by flipping the magnification changing lever from left to right or right to left. The overall magnification power is depended on the selection of the eyepieces (10x or 15x).

Objective lens	1x	1.6x	1x	1.6x
Eyepiece	10x	10x	15x	15x
Total Magnification	10x	16x	15x	24x
Magnification field of view (in mm diameter)	20	12.5	14	9

Preparing the patient



For general hygiene, IBEX recommends the use of a fresh disposable chinrest paper for every patient. Additionally, IBEX recommends cleaning of the head band before each patient exam.

Insure that the illumination intensity knob on the transitional base of the slit lamp is set to the lowest level after each patient exam to prevent uncomfortably bright light reaching the next patient.

Insure the eye pieces have been adjusted in accordance with the examiner.

Patient Alignment

Arrange the slit lamp so your patient will be in a comfortable position. For a quality eye exam, ensure the patient's forehead maintains contact with the forehead rest at all times.

Adjust the chin rest on the headrest assembly up or down by using the chin rest elevation knob so that the patient's eyes are approximately level with the headrest guide mark.



Slowly adjust the intensity knob until a dim light is visible.

Select aperture size 14 then turn the slit width control knob so the light beam is in a thin long vertical line.

Rotate the joystick, until the light beam is centered on the patient's pupil.

Basic Instrument Operation



This coarse alignment should be carried out with the user's naked eye.

Holding the Joystick inclined towards yourself, move/slide the slit lamp transitional base left or right and forward, slowly, closer to the patient, by sliding and rolling the unit on its slide bar (axle) until the slit appears to be approximately in focus on the cornea. Fine adjustment or focus is achieved by carefully tilting the Joystick forward, toward the patient, while observing the slit image through the Optical head.

Maintenance & Cleaning

Routine maintenance



Never attempt to connect or disconnect this instrument with **wet or damp** hands.

Do not disconnect this instrument by pulling on the power cord. Grip and hold the connector firmly when disconnecting from the power source.



Before performing any maintenance or cleaning procedure, ensure the rocker switch is set to the "O" off position.

IBEX LEDs have a life expectancy of 10,000+ hours and should not require changing by the user. We recommend turning off the device after each use to conserve energy and LED life.

In the unlikely event of failure of the LED, please contact your IBEX representative for guidance.

Fuse Replacement

The power supply fuses are located at the rear of the power supply assembly.



To help prevent the possibility of electrical shock ***always*** unplug the slit lamp from its power source before removing or replacing a fuse.



1. First disconnect the slit lamp from the main power source by unplugging the power cord.
2. Insert a small flathead screwdriver into one of the slots and gently pry open the fuse holder located on the back on the power supply.



3. Replace the fuse with the same kind (1Amp 250Volt Slow Blow).
4. Push the fuse holder back in until a “Click” sound is heard and the fuse holder is locked in place.

General cleaning



Before any cleaning of the instrument or the base unit insure that the power cord is disconnected from the power source.

To keep your slit lamp in optimum condition, the slit lamp should be cleaned with a **lightly damp lint free cloth** whenever dust or dirt is detected. Do not use any form of chemical liquids and/or corrosive agents to clean the slit lamp.

Tips for taking care of your optical parts

- Since oil and debris from your hands or from used lens tissue can stain or damage optical coatings, you should not touch any transmissive or reflective surface of your optical system and never reuse a lens tissue.
- *If it's not dirty, don't clean it!* Handling optics increases their chances of getting dirty or damaged, so you should clean optics only when necessary.
- Inspect the eyepieces for dust and stains by holding it near a bright visible-light source. Viewing the optic at different angles allows you to see any scattering of dust and stains.

Recommended Lens Cleaning Technique and Procedure

Step 1: Remove dust

Dusting is always the first step in cleaning your optics. Wiping a dusty optic is like cleaning it with sandpaper. So always dust with a canned air duster, soft optical brush or air blower before wiping any optic. If the dusted optic has no visible stains after you dust it, remember: *"If it's not dirty, don't clean it."* If it's still not clean, proper use of solvents and lens tissue can often do the trick.

Step 2: Wiping the optic

Remember "always wipe slowly and clean the edges first."

Blow off the dust before attempting to clean any optics. Place the brush on the optic surface, apply slight pressure and slowly wipe straight across, from one edge of the optic surface to the other. If the optic still appears dirty proceed to step 3.

Step 3 - Use of solvents and lens tissue

Glass-cleaning solvents will streak, and tissue paper or a t-shirt can scratch, so always clean optics with optical cleaning solution and a low-lint tissue manufactured for cleaning optics. Always use lens tissue with a solvent, because dry lens tissue can scratch optical surfaces. Cleaning your optic's edges before cleaning its faces prevents dirt from being drawn up onto the face. Wiping slowly allows the solvent to evaporate without streaking. Remember, slow and steady cleans the optic.

Technical Specifications

Microscope	Converging- Binocular	Input Voltage	110V/220V AC, 50/60 Hz
Magnification Change Eye Piece	2 Position Objective Lever 10x or 15x	LED Input Filters	3.4V, 700mA Yellow, Cobalt blue, Red-free (Green)
Magnification Ratio	10x, 15x, 16x, 24x	Power Consumption	2.4W (Max)
Real field of view in mm	20, 12.5	Longitudinal (In/Out)	3.9in. (99mm)
Pupillary Distance	49 – 90mm Adjustable	Lateral (Left/Right)	4.6in. (118mm)
Diopter Adjustment	+/- 8D	Vertical (Up/Down)	1.18in. (30mm)
Working Distance	100mm	Chin Rest Range	4.7in. (120mm)
Slit Width	0 – 14mm Continually	Dimension (inch)	28.5(L) x 18.0(W) x 21.9(H)
Slit Length	1 – 14mm Continually	Net Weight	44lb (20kg)
Slit Angles	0° - 180°		

General Classifications

Mode of operation	Intermediate
Degree of mobility	Portable
Type of protection against Electrical Shock	Class I
Degree of protection against Electrical Shock	Type B
Protection against harmful ingress of water	IPX-0
Classification	CE-Regulation 93/42/EEC

Basic Troubleshooting



To help prevent the possibility of electrical shock always unplug the power source.


If the main illumination lamp does not light, check the following:

1. Check the ON/OFF switch located on the front of the unit
2. Turn the intensity knob clockwise for brightness.
3. Check the slit width and aperture to ensure it is opened.



Due to the complicated design of the MCPCB (LED star). Only an authorized Trevi Technology® technician can perform the procedure of replacing the LED star. Please contact your dealer for further instruction.

Block diagram

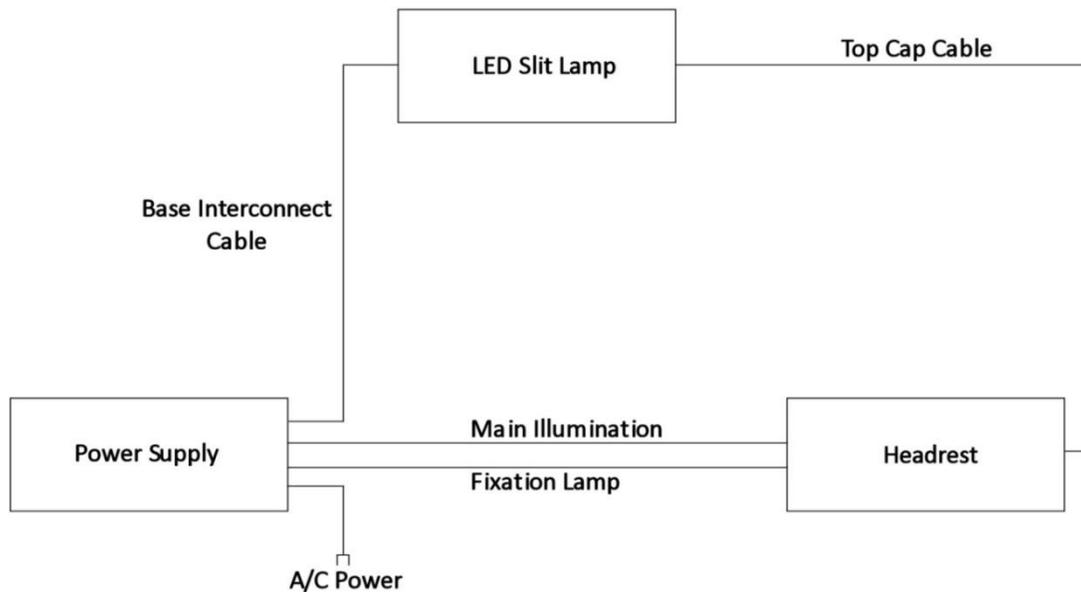

Storage, Transportation & Environment

Storage: Temperature -5°C to +50°C
 Air Pressure 35% to 95%
 Relative Humidity 700kPa to 1060kPa

Transport: Temperature -40°C to +70°C
 Air Pressure 35% to 95%
 Relative Humidity 700kPa to 1060kPa

Environment Condition of Use:

Temperature 0°C to + 50°C
 Air Pressure 35% to 95%
 Relative Humidity 700kPa to 1060kPa

Warranty

Three-Year Comprehensive Warranty

IBEX Slit Lamps are guaranteed for thirty-six months against faulty workmanship, materials or factory defects. The Warranty is a Return-To-Base (RTB) basis at the cost of the customer. All parts and services rendered are included at no additional cost to the customer and in most cases the work is completed within 2-3 business days.

The **IBEX** LED illumination system is designed to operate with consistency for 10,000+ hours of continuous use and therefore is considered as a non-consumable item. In the unlikely event of an LED failure within or after the warranty period, please contact Trevi Technology, Inc. for a replacement.

In case of any defect during warranty period, the concerned person / hospital / institution, shall immediately contact Trevi Technology®. Trevi Technology® will rectify the problem and replace any parts if necessary with no cost.

Quality Systems

Safety Instructions



Read and follow the instructions, cautions and warnings before installing or using this instrument.

The repair/service of this instrument will only be performed by a qualified technician authorized by Trevi Technology®.

Trevi Technology® is responsible for the safety and performance of the equipment only if it is used in accordance with its intended use and the instructions given in the user guide.

This instrument is designed to perform as per the declared indications for safe and reliable service.

Do not modify and/or repair this equipment without authorization of the manufacturer. The manufacturer declines any and all responsibility for loss and/or damages resulting from unauthorized repairs and/or modifications. Any repairs must only be performed by trained and authorized specialists.

Check your slit lamp periodically for any signs of damage or misuse.

This device is intended to be used only by suitably trained and authorized healthcare professionals. Federal law restricts this device to sale by or on the order of a licensed physician.



Never use this device in a potentially explosive environment where volatile solvents (alcohol, petrol, etc.) or flammable anesthetics are in use.



Connect the equipment to properly powered outlets. The entire electrical system in the room or building where your slit lamp is located should conform to IEC and/or local government regulations.

Disconnect the equipment from AC power before servicing and/or cleaning.

Never touch the cable and main power supply system with **wet or damp** hands.

Before any cleaning of the instrument or the base unit insure that the power lead is disconnected.



LED fixtures can reach high temperatures during use – allow ample time to cool before handling any potentially hot surface.

Safety Standards

The FDA and CE marks on the **IBEX** Slit Lamp indicates that it has been tested to meet and/or exceed the U.S. Department of Health and Human Services, Food and Drug Administration Guidance for Slit Lamp Optical Equivalency and Radiation Safety, International Standards within the 93/42/EEC Medical Device Directive and Electrical Safety.



The optical radiation emissions from the slit lamp do not exceed the Threshold Limit Values (TLVs) for optical radiation, as established by the American Conference of Governmental Industrial Hygienists (ACGIH) under worst case clinical exposure conditions and times for ultraviolet and infrared radiation, and visible and near infrared radiation including blue light and aphakic hazards.

It is recommended that the intensity of the light directed into the patient's eye be limited to the minimum level necessary for diagnosis. IBEX Slit Lamps are equipped with protective filters that eliminate UV radiation and short wavelength blue light. Children, infants, aphakes and any person with diseased eyes may be at greater risk. The risk may also be increased if the person being examined has had any exposure with the same instrument or any other ophthalmic instrument using a visible light source during the previous 24 hours. This will apply particularly if the eye has been exposed to retinal photography.

The **IBEX** Slit Lamp is manufactured under strict quality manufacturing systems and certifications as established by ISO 9001:2008 and ISO 13485:2012. These systems encompass materials, processes, assembly and packaging.

Standards Compliance

IBEX Slit Lamp development, production, testing, set-up, maintenance and repair is compatible with IEC 60601-1 (electrical safety), IEC 60601-1-2 (electromagnetic compatibility), EN ISO 15004-2 (optical radiation safety). In addition, **IBEX** Slit Lamps are compliant with the Directive 93/42/EEC as confirmed by the CE marking.

This device meets all the requirements for electromagnetic compatibility according to IEC 60601-1-2. This device is built so that the generation and emission of electromagnetic interference is limited to the extent that other devices are not disturbed in their use in accordance with the regulation and so that the slit lamp itself has appropriate immunity to electromagnetic interference.

This device is designed to comply with the following safety aspects.

Leakage current is less than 0.5mA in normal condition.

Leakage current is less than 1.0mA in a single fault condition.

Isolation voltage between ground and main wires is greater than 1500V.

Contact Information

If for some reason your, **IBEX 2-Step LED Slit Lamp - Compact** does not perform properly or you incur a problem with the instrument:

- a. Consult the applicable section of this manual.
- b. Check: www.ibexeye.com
- c. E-Mail: info@ibexeye.com
Attn: Engineering Department
- d. Call: 614-754-7175
Ask for the Engineering Department



IBEX

Trevi Technology, Inc.
1029 Dublin Rd.
Columbus, OH 43215
614-754-7175
www.ibexeye.com